

X-Ray Protective Headwear



X-Ray Protective

Ongoing Discussion

An increased risk in the occurrence of brain tumors in medical personnel due to high dose work environments, such as interventional radiology, is currently a focus of research.

Particularly in work places where dose values cannot be reduced by a radiation protective shield or another protective system, the head is exposed to a significant amount of radiation.

The figures 1 and 2 exemplify typical medical procedures involving X-Ray radiation without a sufficient radiation protective shield.





The Challenge

To fulfil the high demands for radiation protection of the head and to ensure a comfortable usage, a thoroughly thought out design of the protective headwear is required.

The already required X-Ray protective clothes for medical personnel in high dose work environments create an additional physical exertion. By wearing protective headwear as a supplementary part of the personal protection, the physical stress for the user is increased.

Solutions for Experts

Keeping the needs of our customers in mind, MAVIG developed three individual X-Ray protective headwear solutions. The three designs have different protective and air permeable surface areas.

The design protected variations can be optimally adjusted to individual head circumferences. By this one-size-fits-all approach, the protective headwear features a correct fit and high level of wearing comfort for the user.

Headwear Features

Material Weight-optimized Novalite X-ray protective material

Lead Equivalent Value Pb 0.25 mm, Pb 0.35 mm, Pb 0.50 mm

Outer Cover / Colour ComforTex® HPMF available in the colours:

Curacao, Regatta, Ocean, Lagoon, Limette, Indian Summer, Orchidee

Head Sizes Suitable for head sizes 54 - 63 cm.

Embroidery Embroidery on the left side upon request.



MAVIG X-ray protective headwear complies with IEC 61331-1:2014 / DIN EN 61331-1:2016.

Naturally, these models also have passed EC-type examination and are certified in accordance to PPE 89/686/EEC (CE/PPE).

MAVIG X-ray protective clothing provides the assurance of fulfilling all applicable legal specifications.

Headwear for Medical Personnel

Series RA610 Protective Head Band

The head band protects the front and side of the head. The top of the cranium is not covered by the head band.

Nevertheless, the head band offers an additional X-Ray shielding of the cranium in depending on the angle incident.

Series RA611 Protective Head Cap with Breathable Airflow Strip

This model is the perfect compromise between an enhanced radiation protection of the cranium and heat transmission.

Series RA612 Enclosed Protective Head Cap

The enclosed head cap provides radiation protection of front, side, and top of the head.

Heat transmission is lessened.




















PROTECTIVE EFFECT

The choice of using a 0.25 Pb mm lead equivalent was made after much thought in finding a correct balance between protection and weight

(Please see the comparison chart for more information)

Higher lead equivalent values are of course feasible, but this woul

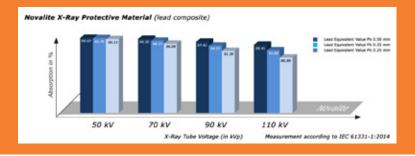



MAVIG GmbHHeadquarters

PO Box 82 03 62 81803 Munich Germany

Stahlgruberring 5 81829 Munich Germany

Phone: +49 (0) 89 420 96 0 Fax: +49 (0) 89 420 96 200 e-Mail: info@mavig.com

Nordic & Baltic Countries

MAVIG Nordic

Stockholm Sweden

Phone: +46 (0) 722 25 25 68 e-Mail: larsson@mavig.com

Benelux, GB, Ireland MAVIG B.V.

Mercuriusweg 86 2516AW Den Haag Netherlands

Phone: +31 (0) 70 33 11 688 Mobile: +31 (0) 61 595 43 48 e-Mail: simmonds@mavig.nl

France MAVIG France SARL

66, Ave. des Champs Elysées F-75008 Paris France

Phone: +33 (0)1 30 59 46 23 Fax: +33 (0)1 30 59 46 23 e-Mail: info@mavig.fr

China MAVIG Healthcare Beijing Co., Ltd.

Block C - Room No. 2607 Wangjing SOHO Tower 1 Futong East Street No. 1 Wangjing Area, Chaoyang District Beijing 100102, China

Phone: +86 - 10 - 64 78 19 21 Fax: +86 - 10 - 64 78 19 21 Mobile: +86 - 13 70 122 85 86 e-Mail: zhang@mavig.com

REPRESENTATION

USA & Canada

Ti-Ba Enterprises, Inc.

25 Hytec Circle Rochester, NY 14606 USA

Phone: +1 (1) 585 247 1212 Fax: +1 (1) 585 247 1395 (international) 1 800 836 8422 (in the USA)

e-Mail: mavigusteam@mavig.com

Middle East & Northern Africa

Mena Medical Development

Ashrafieh – Adlieh Square – Alfaras St. - Alboustany Building, 5th floor Beirut Liban

oan oone: +961 14 23

Phone: +961 14 23 499 Fax: +961 14 26 499 e-Mail: mavig@mena-md.com

www.mavig.com